



ViscoLab PVT+

High-Pressure Viscosity Analysis

- 🌐 Maximizes Instrument Yield with Widest Operational Range and Minimal Sample Volume Requirement
- 🌐 Highest Level of Accuracy and Reliability Enable Key Decision Making
- 🌐 Small Footprint, Turn-key Solution Assures a Straightforward, User-Friendly System Operation
- 🌐 Flexible PVT+ Configurations Offer Unmatched Pressure and Temperature Capabilities

ViscoLab PVT+

FAST AND RELIABLE VISCOSITY MEASUREMENT

For oil, gas, polymers, and supercritical fluids, ViscoLab PVT+ delivers the highest levels of accuracy, speed, and reliability. In high-pressure processes, correctly analysing precious samples is as challenging as it is critical. With millions of dollars on the line, “close enough” is never good enough, and “soon enough” is never fast enough. Since viscosity is an important factor in the quality and marketability of the fluid, measurements have to be accurate and reliable each and every time.



GOLD STANDARD IN HIGH PRESSURE VISCOSITY MEASUREMENT

The Cambridge Viscosity’s ViscoLab PVT+ (Pressure-Volume-Temperature) is the global laboratory standard for determining viscosity at high temperatures and pressures. The ViscoLab PVT+ high-pressure viscometer is designed for viscosity and temperature testing in reservoir fluid analysis, phase behavior, supercritical fluids, and other research and development applications.

Safe and cost-effective to own and operate, the ViscoLab PVT+ requires only 6 ml of sample. Our technology provides statistical certainty that ensures sample conditions are stable, accurate, and repeatable for liquid samples as well as gas and gas condensates.

With its extraordinary capabilities, the ViscoLab PVT+ is the ideal choice for the most critical and challenging applications.

APPLICATION RANGE



- Viscosity measurement from 0.02 to 20,000 centipoise
- Temperatures up to 190°C and pressures up to 30,000 psi (2068 bar) (higher upon request)

STANDARD METHODS



- Complies to:
ASTM D7483
- Correlates to:
ASTM D445


INDUSTRIES



- Oil exploration
- Lube oil development
- Universities
- Research and development
- Chemical manufacturing
- Specialty gases

STREAMLINE REPORTING

Our intuitive HMI integrates time-stamped viscosity, pressure, temperature, density, and shear rate data into a graphical format for easy analysis and storage. All data can be exported over USB. Multi-Language interface available, including Spanish, Portuguese and Chinese. An accessory kit is included with each system to facilitate sample loading, system cleaning, and calibration verification.

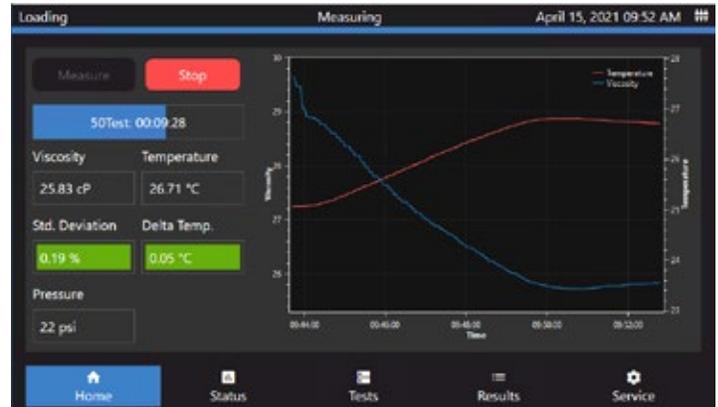

WIDEST OPERATIONAL RANGE IN THE MARKET

The ViscoLab PVT+ features extensive instrument capabilities that enable measurement under challenging downhole conditions. It has the widest operational range in the market:

- Widest pressure range up to 30,000 psi (2068 bar) (higher pressures upon request)
- Widest temperature range up to 190°C (374°F)
- Widest analysis range: 0.02 to 20,000 centipoise (cP)

SUPERIOR SENSOR

Cambridge Viscosity sensors feature an innovative electromagnetic principal that drives an oscillating piston within a precise measurement chamber. Tracking the piston location and travel time gives a direct measure of the test sample viscosity. The constant piston motion generates a continuous data stream which statistically validates the accuracy of the results and presents it in graphical form. There is no need for frequent calibration or maintenance.

SMALL SAMPLE SIZE

Only 6 ml are required in the ViscoLab PVT+ to enable highly accurate measurements for hard to obtain samples whether in reservoir fluid analysis, polymer research, or lube oil/bearing analysis.

HIGHLY ACCURATE AND RELIABLE

- Accurate and reliable measurements on single samples
- Pressure and temperature independently verified
- Analyses liquids and gases
- Optional expansion chamber and adapter for measuring flowing samples

FAST, CONVENIENT, AND LOW COST OF OWNERSHIP

- Small footprint saves laboratory space
- Temperature stability and test completion in 45 minutes or less
- Mercury-free, RoHs compliant and low power consumption

SPECIFICATIONS

Overall viscosity	0.02 To 20,000 centipoise (cP)
Piston range	50:1 Dynamic range
Viscosity accuracy	± 1.0% for liquid viscosities 0.2 to 20,000 cp
Repeatability	± 0.8% of reading
Temperature sensor	PT100
Wetted materials	17-4PH SS or Inconell 718 and C276 Hastelloy
Maximum particle size	25 - 360 Microns
Maximum temperature	190°C/374°F
Maximum pressure	30,000 psi (2068 bar) (Higher upon request)
Power	100 - 240 VAC/50 - 60Hz

PRINCIPLE OF OPERATION

The System At A Glance

Designed for ease of use and accuracy, the complete ViscoLab PVT+ system includes an advanced sensor mounted on a sample conditioning platform with valves, tubing, and pressure transducer. Temperature is controlled from our intuitive interface via an oil bath.

How It Works

The ViscoLab PVT+ combines our next generation electronics with an advanced Cambridge Viscosity SPL-440 sensor. The system employs an integrated recirculating bath that controls temperatures with minimal warm-up time. Sample flow is controlled with a simple three-valve plumbing configuration.

ABOUT CAMBRIDGE VISCOSITY

With more than 10,000 installations worldwide, Cambridge Viscosity is the proven leader in viscosity management technology. With over 25 years of experience, Cambridge Viscosity understands and meets the needs of laboratory researchers and process engineers in a wide range of industries whose jobs depend on the quality, accuracy, and reliability of viscosity measurement equipment. With their patented sensor technology, Cambridge Viscosity has become the gold standard in small sample viscosity measurement.



Contact us for more details.
Visit our website to find the representative closest to you.

ABOUT PAC

PAC develops advanced instrumentation for lab and process applications based on strong **Analytical Expertise** that ensures **Optimal Performance** for our clients. Our analyzers help our clients meet complex industry challenges by providing a low cost of ownership, safe operation, high performance with fast, accurate, and actionable results, high uptime through reliable instrumentation, and compliance with standard methods.

Our solutions are from industry-leading brands: AC Analytical Controls, Advanced Sensors, Alcor, Antek, Herzog, ISL, Cambridge Viscosity, and PetroSpec. We are committed to delivering superior and local customer service worldwide with 16 office locations and a network of over 50 distributors.

HEADQUARTERS

PAC LP | 8824 Fallbrook Drive | Houston, Texas 77064 | USA
T: +1 800.444.8378 | F: +1 281.580.0719



Contact us for more details.
Visit our website to find the PAC representative closest to you.